


**New**  
**Metrowax No.8**  
Ask for details



# Decades of Manufacture

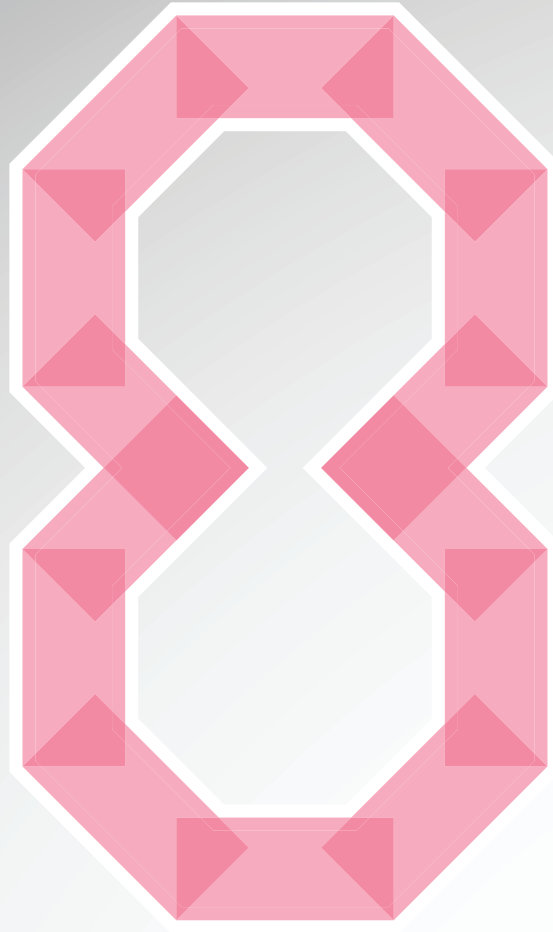

**Since 1936 Metrodent has manufactured quality dental products including dental waxes, denture base and duplicating gels. The Metrodent brand boasts over 80 products in its range, many being market leaders around the world.**



# Content




British Dental Technology  
Official Sponsor

Page 6	Bite Blocks
Page 8	Boxing In Wax
Page 10	Gnathic Waxes
Page 12	Ribbon Wax
Page 14	Soft Bite Wax
Page 16	Strip Wax
Page 18	Wax Up Sticks
Page 22	Metrowax - No.1
Page 24	Metrowax - No.2
Page 26	Metrowax - No.4
Page 28	Metrowax - No.6
Page 30	Metrowax - No.8
Page 32	Metrowax - Light Pink
Page 34	Metrowax - Metromorphic
Page 36	Metrowax - College
Page 40	Metrocryn Denture Base - ST Veined
Page 41	Metrocryn Denture Base - Translucent Veined
Page 42	Metrocryn Denture Base - Standard Pink
Page 43	Metrocryn Denture Base - Translucent Non Veined
Page 44	Metrocryn Denture Base - Natural
Page 45	Metrocryn Denture Base - Ethnic
Page 46	Metrocryn Denture Base - Dark V1 Veined
Page 47	Metrocryn Denture Base - Light Veined
Page 48	Metrocryn Denture Base - Clear
Page 50	Metropair - Denture Repair Acrylic
Page 52	Metrocryn HI (High Impact) Denture Base - Light Pink Veined
Page 54	Orthomet - Orthodontic Acrylic
Page 55	PCS - Plaster Coating Solution
Page 56	Metroform - Custom Impression Tray Acrylic
Page 57	Gypsolve - Cleaning Fluid
Page 58	Metroclens - Cleaning Fluid
Page 59	Duplicating Gels, Metrogel/Supergel/Plastergel
Page 62	Steribim - Germicidal Pumice Compound
Page 64	Metroclave/Hi-Lo - Cold Cure Vessels
Page 65	Flange Mouthpacks

**With over 80 years manufacturing experience  
Metrodent produce a comprehensive range of  
dental waxes tried and trusted by dentists and  
dental technicians both in the UK and around the  
world. Using the finest raw materials Metrodent  
waxes are blended and produced to strict quality  
control standards giving you a first class British  
made product.**



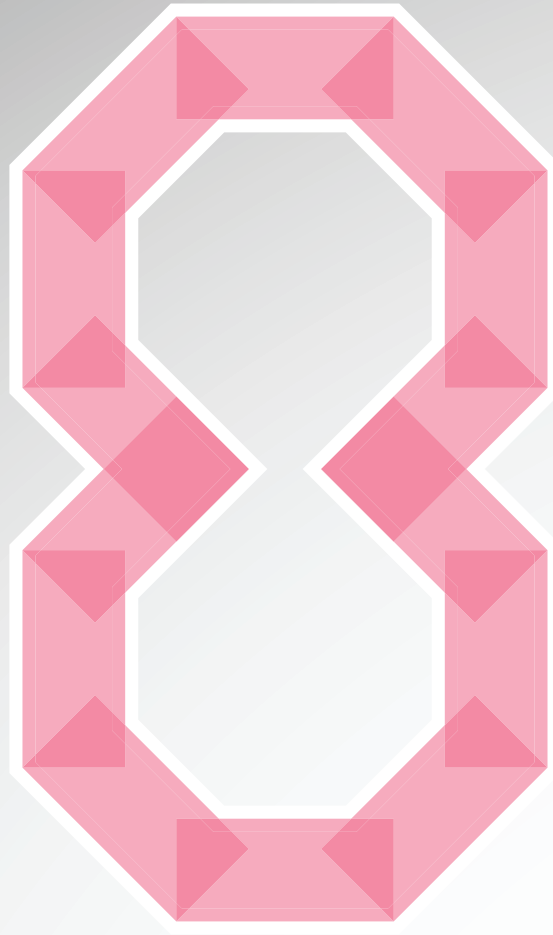

# Types of Dental Waxes



**Made by Metrodent.**

**Bite Blocks**

**Boxing In Wax**

**Gnathic Wax**

**Modelling Wax**

**Ribbon Wax**

**Soft Bite Wax**

**Strip Wax**

**Wax Up Sticks**

# Bite Blocks

## For detailed occlusal records

Metrodent Bite Blocks can be easily adapted to shape and carved with a wax knife. Designed to represent the tooth position and are easy to use. Ideal for a wide range of temperature and climate conditions.

<b>Colours</b>	<b>Pink</b>
	<b>Light Pink</b>
<b>Block Size</b>	<b>114 x 13 x 10 mm</b>
<b>Contents</b>	<b>48 Blocks (8 x 6)</b>
<b>Pack Size</b>	<b>585 g</b>
<b>Pink Order Ref:</b>	<b>WAXBB48</b>
<b>Light Pink Order Ref:</b>	<b>WAXBB48LP</b>



Made by Metrodent





# Boxing In Wax

## For boxing impressions prior to casting

Suitable for production of upper and lower jaws, for opposing models as well as individual split-cast bases. It surrounds the casting, can be easily adapted to its form and builds the edge for the model after the pouring of gypsum.

Colour	Red
Sheet Size	300 x 40 x 2.5 mm
Contents	18 Sheets
Pack Size	450 g
Order Ref:	BIW



Made by Metrodent





# Gnathic Wax

## Differentiate your private work at the bite rim stage Available in Sheet Wax and Bite Blocks

Ivory wax sheets and preformed Ivory wax bite blocks for superior occlusal analysis.  
More accurately represents the dentate area of the rim.  
Improves patient communication and the intra oral contrast  
makes the registration process easier.

### Gnathic Sheet Wax

Colour	Ivory
Dimensions	183 x 87 x 1.75 mm 183 x 87 x 1.50 mm
Pack Size	500 g
Order Ref: 1.75 mm	GWAX500
Order Ref: 1.50 mm	GWAX500EXPORT

### Gnathic Bite Blocks

Colour	Ivory
Dimensions	114 x 13 x 10 mm
Contents	48 Blocks (8 x 6)
Pack Size	585 g
Order Ref:	GWAXBB48









# Ribbon Wax

## For dorsal damming and other intra oral applications

Metrodent Ribbon Wax also known as Periphery Wax is used by dentists and technicians to extend the peripheral edge on an impression tray. High degree of plasticity and exact dorsal dam.

Colour	Red
Contents	126 Strips
Strip Size	210 x 8.0 x 1.50 mm
Pack Size	400g
Order Ref:	RWAX



Made by Metrodent



# Soft Bite Wax

## Economical bite registration wax

Metrodent Soft Bite Wax is a proven bite registration material, formulated to provide dentists with a dense yet easily softened bite wax that is stable at room temperature. The thickness has been specially chosen to provide the most accurate registration.

<b>Colour</b>	<b>Red</b>
<b>Contents</b>	<b>24 Strips (6 x 6 sheets)</b>
<b>Strip Size</b>	<b>230 x 12 x 7.5 mm</b>
<b>Pack Size</b>	<b>400 g</b>
<b>Order Ref:</b>	<b>SRBW</b>



Made by Metrodent





# Strip Wax

## An alternative bite registration wax

Metrodent Strip Wax is a proven bite registration material, formulated to provide dentists with a dense yet easily softened bite wax that is stable at room temperature.

The thickness has been specially chosen to provide the most accurate registration.

<b>Colour</b>	<b>Red</b>
<b>Contents</b>	<b>35 Strips (7 x 5 sheets)</b>
<b>Strip Size</b>	<b>230 x 12 x 7.5 mm</b>
<b>Pack Size</b>	<b>200 g</b>
<b>Order Ref:</b>	<b>SWAX</b>



Made by Metrodent





# Wax Up Sticks

## For setting up and waxing try-in dentures

Metrodent Wax Up Sticks are specially formulated for the waxing up process of either the bite block or try in stages. It carves easily and boils out cleanly.

### Colours

Pink

### Contents

Light Pink

### Strip Size

64 Sticks

### Pack Size

ø 9 x 116 mm

### Order Ref: Pink

400 g

### Order Ref: Light Pink


WUS

LPWUS



**Metrowax stretch toughened dental modelling waxes have been crafted by Metrodent for over 8 decades. Over time production techniques have evolved resulting in products of the highest quality and with the introduction of new blends has given technicians a range of waxes tailored to their needs. Where ever you are, Metrodent has a Metrowax that works for you.**



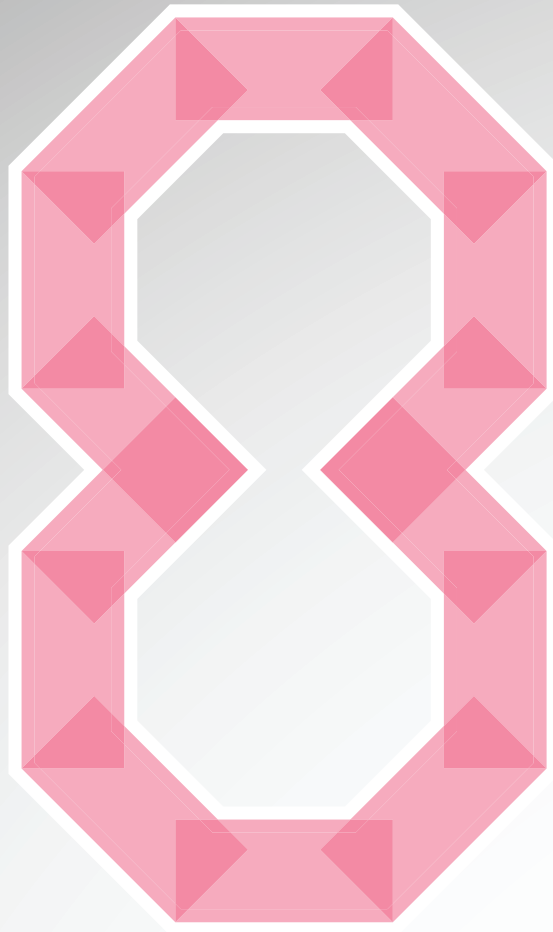

# Modelling Waxes



**Made by Metrodent.**

**Metrowax No.1**

**Metrowax No.2**

**Metrowax No.4**

**Metrowax No.6**

**Metrowax No.8**

**Light Pink**

**Metromorphic**

**College**

# Metrowax No.1

## Highest melting point, ideal for carving

Our highest melting point wax originally designed for warmer climates. Utilising a higher melting point paraffin wax giving a longer molecular chain, it is more stable at higher ambient temperatures and therefore gives excellent carving properties. Ideally suited for high end prosthetics where the try-in is extensively contoured and detailed and where the try-in may be in the mouth for an extended period. Perfect for a hot laboratory in the summer months No.1 will be more brittle in lower temperatures.

<b>Colour</b>	<b>Very Pale Pink</b>
<b>Sheet Sizes</b>	<b>183 x 87 x 1.75 mm</b> <b>183 x 87 x 1.50 mm</b>
<b>Pack Sizes</b>	<b>500 g</b> <b>2.5 Kg</b> <b>20 Kg</b>
<b>Order Ref: 1.75 mm x 500 g</b>	<b>1WAX500</b>
<b>Order Ref: 1.75 mm x 2.5 Kg</b>	<b>1WAX2.5</b>
<b>Order Ref: 1.75 mm x 20 Kg</b>	<b>WAX20</b>
<b>Order Ref: 1.50 mm x 500 g</b>	<b>1WAX500EXPORT</b>
<b>Order Ref: 1.50 mm x 2.5 Kg</b>	<b>1WAX2.5EXPORT</b>
<b>Order Ref: 1.50 mm x 20 Kg</b>	<b>1WAX20EXPORT</b>







# Metrowax No.2

## Our most popular blend

This is our most popular modelling wax both in the UK and abroad. It is a blend of paraffin and micro-crystalline waxes and is designed to give good overall handling properties for everyday laboratory use. Pale Pink in colour No.2 is the longest manufactured blend in our range. Ideal for high workflow prosthetic laboratories where its consistency aids smooth production.

### Colour

Pale 'Bluish' Pink

### Sheet Sizes

183 x 87 x 1.75 mm

183 x 87 x 1.50 mm

183 x 87 x 3.50 mm

### Pack Sizes

500 g

2.5 Kg

20 Kg

Order Ref: 1.75 mm x 500 g

2WAX500

Order Ref: 1.75 mm x 2.5 Kg

2WAX2.5

Order Ref: 1.75 mm x 20 Kg

2WAX20

Order Ref: 1.50 mm x 500 g

2WAX500EXPORT

Order Ref: 1.50 mm x 2.5 Kg

2WAX2.5EXPORT

Order Ref: 1.50 mm x 20 Kg

2WAX20EXPORT

Order Ref: 3.50 mm x 500 g

DTW500

Order Ref: 3.50 mm x 2.5 Kg

DTW2.5

Order Ref: 3.50 mm x 20 Kg

DTW20



Made by Metrodent



# Metrowax No.4

## Darker Pink variant on our most popular blend

This modelling wax is a colour option of our most popular wax. Its handling properties and indications for use are the same as No.2.

<b>Colour</b>	<b>Darker Pink</b>
<b>Sheet Sizes</b>	<b>183 x 87 x 1.75 mm</b> <b>183 x 87 x 1.50 mm</b>
<b>Pack Sizes</b>	<b>500 g</b> <b>2.5 Kg</b> <b>20 Kg</b>
<b>Order Ref: 1.75 mm x 500 g</b>	<b>WAX500</b>
<b>Order Ref: 1.75 mm x 2.5 Kg</b>	<b>4WAX2.5</b>
<b>Order Ref: 1.75 mm x 20 Kg</b>	<b>4WAX20</b>
<b>Order Ref: 1.50 mm x 500 g</b>	<b>4WAX500EXPORT</b>
<b>Order Ref: 1.50 mm x 2.5 Kg</b>	<b>4WAX2.5EXPORT</b>
<b>Order Ref: 1.50 mm x 20 Kg</b>	<b>4WAX20EXPORT</b>



Made by Metrodent





# Metrowax No.6

## Easy to work with enhanced carving properties

Designed to be easy to work with No.6 is a blend of paraffin waxes, micro-crystalline waxes, natural beeswax and resins. The addition of resins enhances the carving properties of the wax making it suitable for contouring and characterisation of the denture base from try-in through to fit. Suitable for prosthetic laboratories with high workflow who may need to characterise their try-ins and require an all round reliable modelling wax.

<b>Colour</b>	<b>Pale Pink</b>
<b>Sheet Sizes</b>	<b>183 x 87 x 1.75 mm</b> <b>183 x 87 x 1.50 mm</b>
<b>Pack Sizes</b>	<b>500 g</b> <b>2.5 Kg</b> <b>20 Kg</b>
<b>Order Ref: 1.75 mm x 500 g</b>	<b>6WAX500</b>
<b>Order Ref: 1.75 mm x 2.5 Kg</b>	<b>6WAX2.5</b>
<b>Order Ref: 1.75 mm x 20 Kg</b>	<b>6WAX20</b>
<b>Order Ref: 1.50 mm x 500 g</b>	<b>6WAX500EXPORT</b>
<b>Order Ref: 1.50 mm x 2.5 Kg</b>	<b>6WAX2.5EXPORT</b>
<b>Order Ref: 1.50 mm x 20 Kg</b>	<b>6WAX20EXPORT</b>



Made by Metrodent



# Metrowax No.8

## Warm Pink wax with superior aesthetics

Our latest edition to the Metrodent range of modelling waxes, No.8 is a complex blend of Paraffin and micro crystalline waxes resulting in a modelling wax with superb handling characteristics. It is Warm Pink in colour with superior aesthetic for both the patient and clinician.

<b>Colour</b>	<b>Warm Pink</b>
<b>Sheet Sizes</b>	<b>183 x 87 x 1.75 mm</b> <b>183 x 87 x 1.50 mm</b>
<b>Pack Sizes</b>	<b>500 g</b> <b>2.5 Kg</b> <b>20 Kg</b>
<b>Order Ref: 1.75 mm x 500 g</b>	<b>8WAX500</b>
<b>Order Ref: 1.75 mm x 2.5 Kg</b>	<b>8WAX2.5</b>
<b>Order Ref: 1.75 mm x 20 Kg</b>	<b>8WAX20</b>
<b>Order Ref: 1.50 mm x 500 g</b>	<b>8WAX500EXPORT</b>
<b>Order Ref: 1.50 mm x 2.5 Kg</b>	<b>8WAX2.5EXPORT</b>
<b>Order Ref: 1.50 mm x 20 Kg</b>	<b>8WAX20EXPORT</b>



Made by Metrodent





# Light Pink

## A very popular shade

The demand for Metrodent to produce a modelling wax in a light pink shade was overwhelming. Rapidly becoming a best seller our Light Pink modelling wax is based on our No.2 formula. A blend of paraffin and micro-crystalline waxes gives familiar ease of handling and consistent qualities aiding smooth work processing in the laboratory. This desirable colour option offers excellent soft tissue representation at try-in stage and is well received by both patients and clinicians.

<b>Colour</b>	<b>Light Pink</b>
<b>Sheet Sizes</b>	<b>183 x 87 x 1.75 mm</b> <b>183 x 87 x 1.50 mm</b>
<b>Pack Sizes</b>	<b>500 g</b> <b>2.5 Kg</b> <b>20 Kg</b>
<b>Order Ref: 1.75 mm x 500 g</b>	<b>LPWAX500</b>
<b>Order Ref: 1.75 mm x 2.5 Kg</b>	<b>LPWAX2.5</b>
<b>Order Ref: 1.75 mm x 20 Kg</b>	<b>LPWAX20</b>
<b>Order Ref: 1.50 mm x 500 g</b>	<b>LPWAX500EXPORT</b>
<b>Order Ref: 1.50 mm x 2.5 Kg</b>	<b>LPWAX2.5EXPORT</b>
<b>Order Ref: 1.50 mm x 20 Kg</b>	<b>LPWAX20EXPORT</b>





# Metromorphic

## A pleasant aroma and easy to work

Metromorphic is a blend of paraffin wax, microcrystalline waxes and natural beeswax which with the omission of resins gives this wax a pleasant aroma. A useful everyday modelling wax.

<b>Colour</b>	<b>Pale 'Bluish' Pink</b>
<b>Sheet Sizes</b>	<b>183 x 87 x 1.75 mm</b> <b>183 x 87 x 1.50 mm</b>
<b>Pack Sizes</b>	<b>500 g</b> <b>2.5 Kg</b> <b>20 Kg</b>
<b>Order Ref: 1.75 mm x 500 g</b>	<b>7WAX500</b>
<b>Order Ref: 1.75 mm x 2.5 Kg</b>	<b>7WAX2.5</b>
<b>Order Ref: 1.75 mm x 20 Kg</b>	<b>7WAX20</b>
<b>Order Ref: 1.50 mm x 500 g</b>	<b>7WAX500EXPORT</b>
<b>Order Ref: 1.50 mm x 2.5 Kg</b>	<b>7WAX2.5EXPORT</b>
<b>Order Ref: 1.50 mm x 20 Kg</b>	<b>7WAX20EXPORT</b>



Made by Metrodent





# College

## For training and education purposes


Ideal for Dental Schools, Colleges and Universities this wax is made of secondary blends of our primary waxes and varies in colour and handling properties.

<b>Colour</b>	<b>Pink</b>
<b>Sheet Sizes</b>	<b>183 x 87 x 1.75 mm</b> <b>183 x 87 x 1.50 mm</b>
<b>Pack Sizes</b>	<b>500 g</b> <b>20 Kg</b>
<b>Order Ref: 1.75 mm/500 g</b>	<b>CWAX</b>
<b>Order Ref: 1.75 mm/20 Kg</b>	<b>CWAX20</b>
<b>Order Ref: 1.50 mm/500 g</b>	<b>CWAX500EXPORT</b>
<b>Order Ref: 1.50 mm/20 Kg</b>	<b>CWAX20EXPORT</b>



Made by Metrodent

**Over decades Metrodent have developed a range of denture base acrylics for both our UK and Export markets. Metrocyl Denture Base is available in 8 shade options including ethnic and clear. In addition we have an orthodontic resin, Orthomet, for the “spray” technique available in clear and tinted pink. All our heat cure resins conform to Standard BS EN ISO20795-1: 2008**

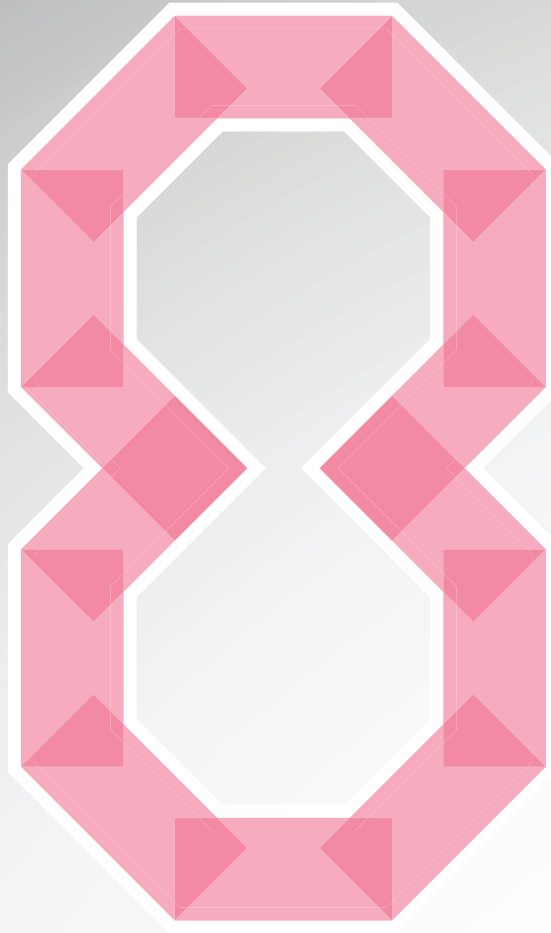

# Shades of Denture Base





**Made by Metrodent.**

**ST Veined**

**Translucent Veined**

**Standard Pink**

**Translucent Non Veined**

**Dark V1 Veined**


**Light Veined**

**Natural**

**Ethnic**

**Clear**

# ST Veined

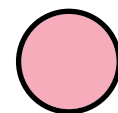

## Pink Translucent Veined

Our most popular acrylic formulated in response to market demand to match the leading brand product. A silky dough is produced with exceptionally good flow properties resulting in reduced flash thickness.

### Universal

- ✓ Snap stage 15 minutes
- ✓ Work time 60 minutes at 23°C
- ✓ Can be quick cured or processed overnight

### Rapid Cure

- ✓ Snap stage 30 minutes
- ✓ Work time 30 minutes at 23°C
- ✓ Can be quick cured in 20 minutes


Order Ref: U1SV	Metrocyl - Universal ST Veined 1 Kg
Order Ref: U3SV	Metrocyl - Universal ST Veined 3 Kg
Order Ref: U12SV	Metrocyl - Universal ST Veined 12 Kg
Order Ref: R1SV	Metrocyl - Rapid Cure ST Veined 1 Kg
Order Ref: R3SV	Metrocyl - Rapid Cure ST Veined 3 Kg
Order Ref: R12SV	Metrocyl - Rapid Cure ST Veined 12 Kg
Order Ref: LU500	Metrocyl - Universal Liquid 500 ml
Order Ref: LU1	Metrocyl - Universal Liquid 1 Litre

Order Ref: MP1SV	Metropair - Universal ST Veined 1 Kg
Order Ref: MP3SV	Metropair - Universal ST Veined 3 Kg
Order Ref: MP12SV	Metropair - Universal ST Veined 12 Kg
Order Ref: MPL500	Metropair - Universal Liquid 500 ml
Order Ref: MPL1	Metropair - Universal Liquid 1 Litre





# Translucent Veined

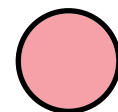

## Medium Pink Veined

A popular translucent denture base with veining.  
Slightly whiter or more opaque than either our Translucent or Light Veined acrylic.

### Universal

- ✓ Snap stage 15 minutes
- ✓ Work time 60 minutes at 23°C
- ✓ Can be quick cured or processed overnight

### Rapid Cure

- ✓ Snap stage 30 minutes
- ✓ Work time 30 minutes at 23°C
- ✓ Can be quick cured in 20 minutes

Order Ref: U1TV	Metrocryl - Universal Translucent Veined 1 Kg
Order Ref: U3TV	Metrocryl - Universal Translucent Veined 3 Kg
Order Ref: U12TV	Metrocryl - Universal Translucent Veined 12 Kg
Order Ref: R1TV	Metrocryl - Rapid Cure Translucent Veined 1 Kg
Order Ref: R3TV	Metrocryl - Rapid Cure Translucent Veined 3 Kg
Order Ref: R12TV	Metrocryl - Rapid Cure Translucent Veined 12 Kg
Order Ref: LU500	Metrocryl - Universal Liquid 500 ml
Order Ref: LU1	Metrocryl - Universal Liquid 1 Litre


Order Ref: MP1TV	Metopair - Translucent Veined 1 Kg
Order Ref: MP3TV	Metopair - Translucent Veined 3 Kg
Order Ref: MP12TV	Metopair - Translucent Veined 12 Kg
Order Ref: MPL500	Metopair - Liquid 500 ml
Order Ref: MPL1	Metopair - Liquid 1 Litre



Made by Metrodent



# Standard Pink

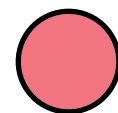

## Natural Pink

With greater opacity this acrylic has no veining.  
A lustrous natural pink mimicking healthy gingival mucosa, un-veined.

### Universal

- ✓ Snap stage 15 minutes
- ✓ Work time 60 minutes at 23°C
- ✓ Can be quick cured or processed overnight

### Rapid Cure

- ✓ Snap stage 30 minutes
- ✓ Work time 30 minutes at 23°C
- ✓ Can be quick cured in 20 minutes

Order Ref: U1S	Metrocyl - Universal Standard Pink 1 Kg
Order Ref: U3S	Metrocyl - Universal Standard Pink 3 Kg
Order Ref: U12S	Metrocyl - Universal Standard Pink 12 Kg
Order Ref: R1S	Metrocyl - Rapid Cure Standard Pink 1 Kg
Order Ref: R3S	Metrocyl - Rapid Cure Standard Pink 3 Kg
Order Ref: R12S	Metrocyl - Rapid Cure Standard Pink 12 Kg
Order Ref: LU500	Metrocyl - Universal Liquid 500 ml
Order Ref: LU1	Metrocyl - Universal Liquid 1 Litre


Order Ref: MP1S	Metopair - Standard Pink 1 Kg
Order Ref: MP3S	Metopair - Standard Pink 3 Kg
Order Ref: MP12S	Metopair - Standard Pink 12 Kg
Order Ref: MPL500	Metopair - Liquid 500 ml
Order Ref: MPL1	Metopair - Liquid 1 Litre



Made by Metrodent



# Translucent Non Veined

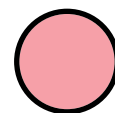

## Light Gingival Pink

Translucent denture base, translucency allows the patient's own gingival tones to "colour" the denture base.

### Universal

- ✓ Snap stage 15 minutes
- ✓ Work time 60 minutes at 23°C
- ✓ Can be quick cured or processed overnight

### Rapid Cure

- ✓ Snap stage 30 minutes
- ✓ Work time 30 minutes at 23°C
- ✓ Can be quick cured in 20 minutes

Order Ref: U1T	Metrocyl - Universal Translucent Non Veined 1 Kg
Order Ref: U3T	Metrocyl - Universal Translucent Non Veined 3 Kg
Order Ref: U12T	Metrocyl - Universal Translucent Non Veined 12 Kg
Order Ref: R1T	Metrocyl - Rapid Cure Translucent Non Veined 1 Kg
Order Ref: R3T	Metrocyl - Rapid Cure Translucent Non Veined 3 Kg
Order Ref: R12T	Metrocyl - Rapid Cure Translucent Non Veined 12 Kg
Order Ref: LU500	Metrocyl - Universal Liquid 500 ml
Order Ref: LU1	Metrocyl - Universal Liquid 1 Litre


Order Ref: MP1T	Metropair - Translucent Non Veined 1 Kg
Order Ref: MP3T	Metropair - Translucent Non Veined 3 Kg
Order Ref: MP12T	Metropair - Translucent Non Veined 12 Kg
Order Ref: MPL500	Metropair - Liquid 500 ml
Order Ref: MPL1	Metropair - Liquid 1 Litre



Made by Metrodent



# Natural

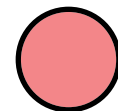

## Medium Pink

A pink opaque denture base mimicking certain gingival tissues. Medium pink with a very subtle orange tint and good opacity.

### Universal

- ✓ Snap stage 15 minutes
- ✓ Work time 60 minutes at 23°C
- ✓ Can be quick cured or processed overnight

### Rapid Cure

- ✓ Snap stage 30 minutes
- ✓ Work time 30 minutes at 23°C
- ✓ Can be quick cured in 20 minutes

Order Ref: U1N	Metrocyl - Universal Natural 1 Kg
Order Ref: U3N	Metrocyl - Universal Natural 3 Kg
Order Ref: U12N	Metrocyl - Universal Natural 12 Kg
Order Ref: R1N	Metrocyl - Rapid Cure Natural 1 Kg
Order Ref: R3N	Metrocyl - Rapid Cure Natural 3 Kg
Order Ref: R12N	Metrocyl - Rapid Cure Natural 12 Kg
Order Ref: LU500	Metrocyl - Universal Liquid 500 ml
Order Ref: LU1	Metrocyl - Universal Liquid 1 Litre


Order Ref: MP1N	Metropair - Natural 1 Kg
Order Ref: MP3N	Metropair - Natural 3 Kg
Order Ref: MP12N	Metropair - Natural 12 Kg
Order Ref: MPL500	Metropair - Liquid 500 ml
Order Ref: MPL1	Metropair - Liquid 1 Litre



Made by Metrodent



# Ethnic

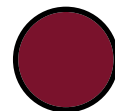

## Dark Purple shade of Pink veined

Developed to match the gingival tissues of patients with afro-caribbean heritage.

### Universal

- ✓ Snap stage 15 minutes
- ✓ Work time 60 minutes at 23°C
- ✓ Can be quick cured or processed overnight

### Rapid Cure

- ✓ Snap stage 30 minutes
- ✓ Work time 30 minutes at 23°C
- ✓ Can be quick cured in 20 minutes

Order Ref: U1E	Metrocryl - Universal Ethnic 1 Kg
Order Ref: U3E	Metrocryl - Universal Ethnic 3 Kg
Order Ref: R1E	Metrocryl - Rapid Cure Ethnic 1 Kg
Order Ref: R3E	Metrocryl - Rapid Cure Ethnic 3 Kg
Order Ref: LU500	Metrocryl - Universal Liquid 500 ml
Order Ref: LU1	Metrocryl - Universal Liquid 1 Litre

Order Ref: MP1EV	Metropair - Ethnic 1 Kg
Order Ref: MP3EV	Metropair - Ethnic 3 Kg
Order Ref: MPL500	Metropair - Liquid 500 ml
Order Ref: MPL1	Metropair - Liquid 1 Litre




Made by Metrodent





# Dark V1 Veined

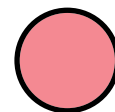

## A more reddish Pink veined

Developed originally for export to our East Asian markets and offering a variation in gingival tissue matching.

### Universal

- ✓ Snap stage 15 minutes
- ✓ Work time 60 minutes at 23°C
- ✓ Can be quick cured or processed overnight

### Rapid Cure

- ✓ Snap stage 30 minutes
- ✓ Work time 30 minutes at 23°C
- ✓ Can be quick cured in 20 minutes

Order Ref: U1DV	Metrocyl - Universal Dark V1 Veined 1 Kg
Order Ref: U3DV	Metrocyl - Universal Dark V1 Veined 3 Kg
Order Ref: R1DV	Metrocyl - Rapid Cure Dark V1 Veined 1 Kg
Order Ref: R3DV	Metrocyl - Rapid Cure Dark V1 Veined 3 Kg
Order Ref: LU500	Metrocyl - Universal Liquid 500 ml
Order Ref: LU1	Metrocyl - Universal Liquid 1 Litre

Order Ref: MP1DV	Metropair - Dark V1 Veined 1 Kg
Order Ref: MP3DV	Metropair - Dark V1 Veined 3 Kg
Order Ref: MPL500	Metropair - Liquid 500 ml
Order Ref: MPL1	Metropair - Liquid 1 Litre



Made by Metrodent



# Light Veined

## Light Translucent Pink Veined

Similar to our most popular ST Veined acrylic but in a lighter more subtle pink suitable for patients with paler gingival tissues.

### Universal

- ✓ Snap stage 15 minutes
- ✓ Work time 60 minutes at 23°C
- ✓ Can be quick cured or processed overnight

### Rapid Cure

- ✓ Snap stage 30 minutes
- ✓ Work time 30 minutes at 23°C
- ✓ Can be quick cured in 20 minutes

Order Ref: U1LV	Metrocyl - Universal Light Veined 1 Kg
Order Ref: U3LV	Metrocyl - Universal Light Veined 3 Kg
Order Ref: U12LV	Metrocyl - Universal Light Veined 12 Kg
Order Ref: R1LV	Metrocyl - Rapid Cure Light Veined 1 Kg
Order Ref: R3LV	Metrocyl - Rapid Cure Light Veined 3 Kg
Order Ref: R12LV	Metrocyl - Rapid Cure Light Veined 12 Kg
Order Ref: LU500	Metrocyl - Universal Liquid 500 ml
Order Ref: LU1	Metrocyl - Universal Liquid 1 Litre


Order Ref: MP1LV	Metopair - Translucent Non Veined 1 Kg
Order Ref: MP3LV	Metopair - Translucent Non Veined 3 Kg
Order Ref: MP12LV	Metopair - Translucent Non Veined 12 Kg
Order Ref: MPL500	Metopair - Liquid 500 ml
Order Ref: MPL1	Metopair - Liquid 1 Litre



Made by Metrodent



# Clear

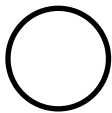

## Clear

Ideal for clear palates, TMJ appliances, orthodontic appliances, surgical guides etc. Using rapid cure gives an exceptionally clear acrylic whereas universal can sometimes have a slight opacity.

### Universal

- ✓ Snap stage 15 minutes
- ✓ Work time 60 minutes at 23°C
- ✓ Can be quick cured or processed overnight

### Rapid Cure

- ✓ Snap stage 30 minutes
- ✓ Work time 30 minutes at 23°C
- ✓ Can be quick cured in 20 minutes
- ✓ Use of rapid cure gives a very clear acrylic\*

Order Ref: U1C	Metrocryn - Universal Clear 1 Kg
Order Ref: U3C	Metrocryn - Universal Clear 3 Kg
Order Ref: U12C	Metrocryn - Universal Clear 12 Kg
Order Ref: R1C	Metrocryn - Rapid Cure Clear 1 Kg
Order Ref: R3C	Metrocryn - Rapid Cure Clear 3 Kg
Order Ref: R12C	Metrocryn - Rapid Cure Clear 12 Kg
Order Ref: LU500	Metrocryn - Universal Liquid 500 ml
Order Ref: LU1	Metrocryn - Universal Liquid 1 Litre

Order Ref: MP1C	Metopair - Clear 1 Kg
Order Ref: MP3C	Metopair - Clear 3 Kg
Order Ref: MP12C	Metopair - Clear 2 Kg
Order Ref: MPL500	Metopair - Liquid 500 ml
Order Ref: MPL1	Metopair - Liquid 1 Litre

\*The use of Universal liquid with clear acrylic will result in yellowing of the cured acrylic and we recommend always using our specially formulated clear liquid.



Made by Metrodent



A close-up photograph of human teeth, showing the upper and lower arches. The teeth are white with visible pinkish-red gums. A semi-transparent pinkish-red overlay is applied to the right side of the image, creating a soft, ethereal effect. The text is overlaid on the left side of the image.

**Made by Metrodent.**  
**ST Veined**  
**Translucent Veined**  
**Standard Pink**  
**Translucent Non Veined**  
**Dark V1 Veined**  
**Light Veined**  
**Natural**  
**Ethnic**  
**Clear**



# Metropair

## Denture Repair Acrylic

Cold cure acrylic for denture repairs, additions, relines etc., non slumping, quick setting and free of porosity if cured in a pressure flask. **Use only with Metropair Liquid.**

Order Ref: MP1SV	Metropair - Universal ST Veined 1 Kg
Order Ref: MP3SV	Metropair - Universal ST Veined 3 Kg
Order Ref: MP12SV	Metropair - Universal ST Veined 12 Kg
Order Ref: MP1TV	Metropair - Translucent Veined 1 Kg
Order Ref: MP3TV	Metropair - Translucent Veined 3 Kg
Order Ref: MP12TV	Metropair - Translucent Veined 12 Kg
Order Ref: MP1S	Metropair - Standard Pink 1 Kg
Order Ref: MP3S	Metropair - Standard Pink 3 Kg
Order Ref: MP12S	Metropair - Standard Pink 12 Kg
Order Ref: MP1T	Metropair - Translucent Non Veined 1 Kg
Order Ref: MP3T	Metropair - Translucent Non Veined 3 Kg
Order Ref: MP12T	Metropair - Translucent Non Veined 12 Kg
Order Ref: MP1DV	Metropair - Dark V1 Veined 1 Kg
Order Ref: MP3DV	Metropair - Dark V1 Veined 3 Kg
Order Ref: MP1LV	Metropair - Light Veined 1 Kg
Order Ref: MP3LV	Metropair - Light Veined 3 Kg
Order Ref: MP12LV	Metropair - Light Veined 12 Kg
Order Ref: MP1N	Metropair - Natural 1 Kg
Order Ref: MP3N	Metropair - Natural 3 Kg
Order Ref: MP12N	Metropair - Natural 12 Kg
Order Ref: MP1EV	Metropair - Ethnic 1 Kg
Order Ref: MP3EV	Metropair - Ethnic 3 Kg
Order Ref: MP1C	Metropair - Clear 1 Kg
Order Ref: MP3C	Metropair - Clear 3 Kg
Order Ref: MP12C	Metropair - Clear 2 Kg
Order Ref: MPL500	Metropair - Liquid 500 ml
Order Ref: MPL1	Metropair - Liquid 1 Litre



Made by Metrodent



# Metrocryn Denture Base



# HI (High Impact)

## Light Pink Veined

Impact strength is different from tensile or flexural strength and will prevent breakage of the denture base due to impact shocks for example if the denture is dropped on a hard surface like a porcelain sink. High impact acrylic is not a substitute for a CoCr denture base where a patient's bite or biting forces would indicate this prescription. In thin cross sections, such as the replacement of a single tooth such as a central or lateral incisor in a close bite situation, high impact acrylics are also no more suitable than a conventional PMMA denture base. Repairs – Should the need arise to repair a denture fabricated from high impact denture base it should be re-processed using heat cure high impact resin to maintain its high impact properties. The use of a standard "cold cure" resin will mean the denture is no longer considered to have greater impact resistance. Unlike heat cure denture base, the high impact powder does not absorb monomer very well due to it being rubber modified. To this end the powder should be weighed and the volume of liquid measured according to the instructions. The powder will need to be worked into the liquid using a metal spatula to ensure the correct amount of powder is used in the mix. Use of the 'pepper and salt' technique is likely to give incorrect polymerisation.

The impact strength of notched and surface defected specimens (according to ISO 179) gives values of 2.55 Kj/m<sup>2</sup> and 9.55 Kj/m<sup>2</sup> respectively.

### Polymer

- ✓ Snap stage 15 minutes
- ✓ Work time 20 minutes at 23°C

Order Ref: HI1

Order Ref: HI3

Order Ref: HI12

Order Ref: HIL500

Order Ref: HIL1

Metrocyl HI - Polymer 1 Kg

Metrocyl HI - Polymer 3 Kg

Metrocyl HI - Polymer 12 Kg

Metrocyl HI - Monomer 500 ml

Metrocyl HI - Monomer 1 Litre

### Monomer

- ✓ Only use with Metrocyl HI




# Orthodontic Acrylic





# Orthomet

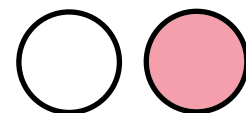

## Clear or Tinted Pink

Orthomet is used in the “spray” technique for the manufacture of Orthodontic appliances. With the use of a pressure flask such as the Metroclave or Hi-Lo, the material will be cured in approximately 10 minutes showing no porosity or shrinkage.

Order Ref: OP1	Orthomet - Polymer 1 Kg
Order Ref: OP3	Orthomet - Polymer 3 Kg
Order Ref: OP10	Orthomet - Polymer 10 Kg
Order Ref: OLC500	Orthomet - Clear Monomer 500 ml
Order Ref: OLC1	Orthomet - Clear Monomer 1 Litre
Order Ref: OLT500	Orthomet - Tinted Pink Monomer 500 ml
Order Ref: OLT1	Orthomet - Tinted Pink Monomer 1 Litre



Made by Metrodent



# Plaster Coating Solution

## Cold Mould Seal

Premium grade formaldehyde free alginate solution for insulation of denture base resins from gypsum moulds and casts when making dentures, acrylic crowns etc.

Order Ref: PCS1	PCS (Plaster Coating Solution) - 1 Litre
Order Ref: PCS2	PCS (Plaster Coating Solution) - 2 Litres
Order Ref: PCS5	PCS (Plaster Coating Solution) - 5 Litres
Order Ref: PCS25	PCS (Plaster Coating Solution) - 25 Litres
Order Ref: PCS200	PCS (Plaster Coating Solution) - 200 Litres



Made by Metrodent

# Metroform

## Custom impression tray acrylic

A high quality self cure material with excellent working properties and economic in use.

Blue - standard time.  
Dough time 2 to 3 minutes.

Green - long work time.  
Dough time 4 to 5 minutes.

Order Ref: MTFB1	Metroform Blue - Polymer 1 Kg
Order Ref: MTFB3	Metroform Blue - Polymer 3 Kg
Order Ref: MTFB12	Metroform Blue - Polymer 12 Kg
Order Ref: MTFG1	Metroform Green - Polymer 1 Kg
Order Ref: MTFG3	Metroform Green - Polymer 3 Kg
Order Ref: MTFG12	Metroform Green - Polymer 12 Kg

Order Ref: MTFL500	Metroform - Monomer 500 ml
Order Ref: MTFL2	Metroform - Monomer 2 Litres



Made by Metrodent



# Gypsolve

## Cleaning Fluid

Gypsolve is a clear non toxic, non corrosive liquid which can safely be used for the rapid removal of residual plaster. Also suitable for use in ultrasonic cleaners.

Order Ref: GG5	Gypsolve - Green 5 Litres
Order Ref: GG25	Gypsolve - Green 25 Litres
Order Ref: GC5	Gypsolve - Clear 5 Litres
Order Ref: GC25	Gypsolve - Clear 25 Litres



Made by Metrodent



# Metroclens

## Cleaning Fluid

Metroclens is a tray cleaner, plaster remover and hard surface cleaner. Used in varying diluted forms gives Metroclens solution optimum cleaning properties. Can be used in an ultrasonic cleaner to speed up the cleaning process. Must not be used for aluminium or aluminium alloys.

Order Ref: MCLENS1

Metroclens - 1 Litre



Made by Metrodent

# Duplicating Gels

## Metrogel

## Supergel

## Plastergel

Metrogel - successful aid in the partial denture technique.

Supergel - for the cobalt chrome partial denture laboratory.

Plastergel - is ideal for the production of quality duplicate models made from plasters and dental stones.

Order Ref: MGEL

Order Ref: MGEL/C

Order Ref: SGEL

Order Ref: SGEL/C

Order Ref: PGEL

Order Ref: PGEL/C

Metrogel - Blue 5 Kg

Metrogel - Clear 5 Kg

Supergel - Green 5 Kg

Supergel - Clear 5 Kg

Plastergel - Red 5 Kg

Plastergel - Clear 5 Kg



Made by Metrodent





**etrogel**

Reversible Agar Hydrocolloid  
Duplicating Gel

Manufactured under licence from BEGO  
Guaranteed to conform to the requirements of EN ISO 14356 2003 type I

Store in a cool dry place

Net Wt.

Minimum volume 4.2 Litres

LOT  
1 D 2017  
2022/04



Metrodent Limited, Lowergate Works, Lowergate, Paddock, Huddersfield HD3 4EP England  
Tel: +44 (0) 1484 461616 Fax: +44 (0) 1484 461700  
sales@metrodent.com metrodent.com

**Supergel**

Reversible Agar Hydrocolloid  
Duplicating Gel

Resilient, superior quality for accurate reproduction of fine detail

Manufactured under licence from BEGO  
Guaranteed to conform to the requirements of EN ISO 14356 2003 type I

Store in a cool dry place

5 Kg Nett Wt.

Minimum volume 4.2 Litres

LOT  
4 C 2017  
2022/02



Metrodent Limited, Lowergate Works, Lowergate, Paddock, Huddersfield HD3 4EP England  
Tel: +44 (0) 1484 461616 Fax: +44 (0) 1484 461700  
sales@metrodent.com metrodent.com

**Plastergel**

Duplicating Gel

Reversible Agar Hydrocolloid  
For plaster and investments. Alum wash not required

Manufactured under licence from BEGO  
Guaranteed to conform to the requirements of EN ISO 14356 2003 type I

Store in a cool dry place

5 Kg Nett Wt.

Minimum volume 4.2 Litres

LOT  
3 C 2017  
2022/02



Metrodent Limited, Lowergate Works, Lowergate, Paddock, Huddersfield HD3 4EP England  
Tel: +44 (0) 1484 461616 Fax: +44 (0) 1484 461700  
sales@metrodent.com metrodent.com







# Steribim

## Germicidal pumice compound

Destroys all germs and the germicidal effect lasts down to the last gram. Special additives make it kind to the skin and prevent the pumice from going sour in the polishing trough. Dust free.

Order Ref: STB12

Steribim Germicidal Pumice - 12 Kg



Made by Metrodent

# Cold Cure Vessels



# Metroclave and Hi-Lo

## The simple and reliable solution

Metrodent have been making the Metroclave and Hi-Lo Cold Cure Vessels for over 30 years. We think they are the best, and many Technicians agree with us. Part of their success is based on their simplicity. To process cold cure resins two things are required heat and pressure. Using the Metroclave and Hi-Lo the heat is supplied by filling the vessels with warm water and the pressure is applied by screwing down the lid and indicated on the gauge. That's it!

### Advantages

- ✓ No compressor required
- ✓ No heating plate needed
- ✓ Does not tie-up the flask press whilst in use
- ✓ Requires little maintenance
- ✓ Gives years of service
- ✓ Inexpensive

**Metroclave:** The lid is recessed in a well in the body of the vessel, this allows it to be used on a bench away from the sink or drainer. If the Metroclave is filled to just the right level, the well allows the lid to be fitted without spillage. Curing time 6 to 15 minutes according to the water temperature. Holds up to 5 models. Dimensions: ø 160 x 90 mm.

**Hi-Lo-** This is smaller of the two vessels and has a capacity for up to three models. It is ideal for the smaller laboratory or department or even for occasional use in a dental surgery. Curing time 6 to 15 minutes according to the water temperature. Dimensions; ø 110 x 110 mm.

**Order Ref: MTC**

**Metrotroclave**

**Order Ref: HLA**

**Hi-Lo**



Made by Metrodent

# Flange Mouthpacks

## Non Sterile

Absorbent sponges to pack the airway in order to prevent inhalation/ ingestion of debris, etc.

Autoclave before use.

Colour	Yellow
Contents	12
Order Ref:	MPACK



Made by Metrodent





